

Compact Embedding Center ECC

Space-saving Paraffin Embedding Center.

Consisting of:

Paraffin dispenser unit with a 4.5 liter paraffin reserve.

Heated tray for cassettes or for entire baskets from almost all manufacturers.

Heated compartment for the pouring molds.

Cooling plate with space for 72 cassettes


Cold spot.

Heated work surface for approx. 20 cassettes.

Drawer to collect excess paraffin.

Adjustable paraffin spout, lighting and optional magnifying glass and / or foot switch.

Touchscreen for setting and monitoring all temperatures, such as programming the operating times.

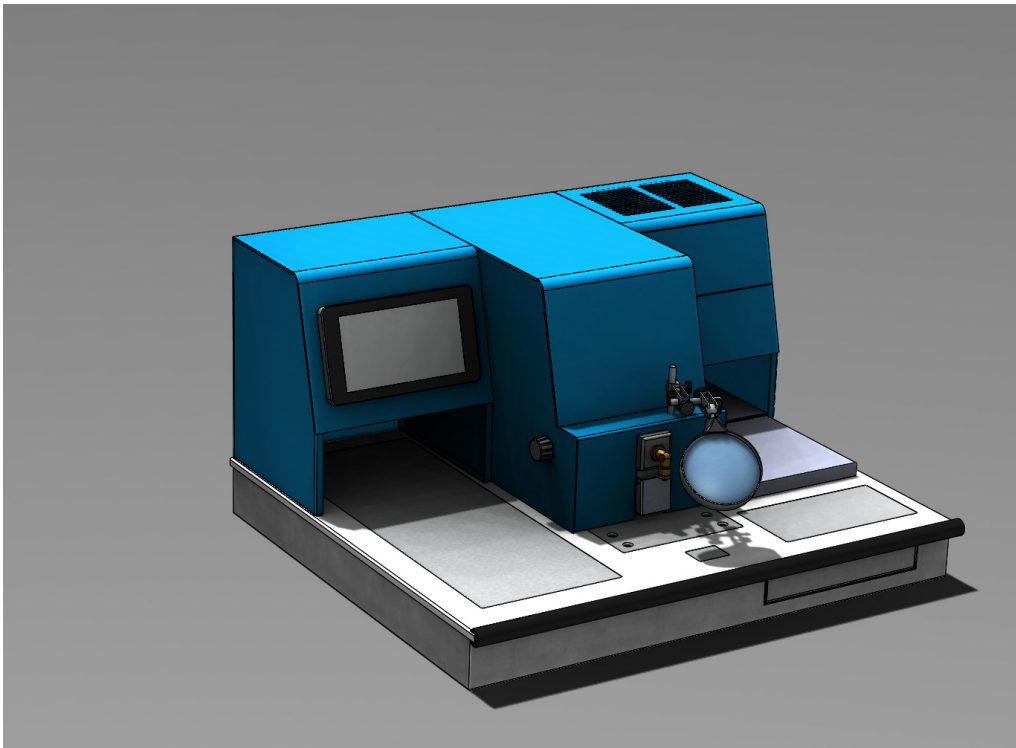

Technical specifications:

Dispenser unit: W630 x D550 x H400mm,

230V / 50Hz / 1KW

Max. 70°C,

4.5 liter paraffin storage tank

Article no .: RWECC