


# Formalin Filling Machine

Formalin Filling Machine for the operation room.

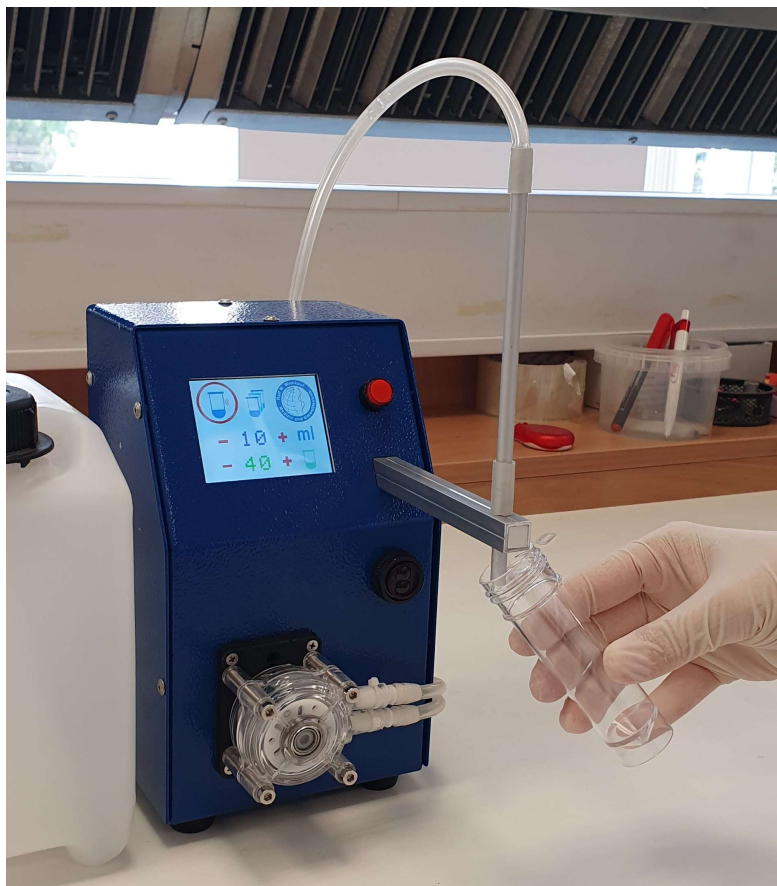

**Non-contact filling of formalin into vessels.**

**The required quantity can be preset on the touch screen, a sensor recognizes the vessel and begins to fill it with the appropriate with the corresponding quantity.**

**The formalin is sucked out of a 5 liter (10 / 20 liter) container.**

**Settings: from 1ml to 10ml in 1ml steps  
from 10ml to 95ml in 5ml steps**

**Accuracy: +/- 2%**

## **Technical data:**

**Dimensions of the device (W/H/D) 210 x 125 x 200mm  
weight: 1,6 Kg  
Power supply: 110V / 60HZ 12V DC power supply unit**

**Article No. 502200**