



## Fast Freezing Cryofix -45°C

### Snap Freezing


Device for rapid freezing of tissue samples.

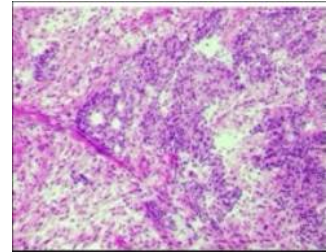

**Liquid Nitrogen**

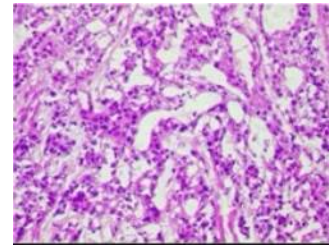

**Dry Ice**

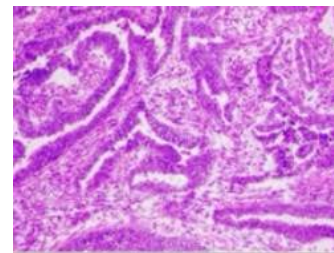

**Isopentane**

## Fast freezing unit Cryofix -45 °C

Instrument for fast freezing of tissue samples.

Fits all conventional cryostats.

Isopentane is cooled in the retort down to -45 °C

The formation of ice crystals almost becomes completely avoided.

Immediate shock freezing reduces the freezing time to a minimum.

No danger as with liquid nitrogen.

The compact design allows installation close to your Cryostat.

Also as an alternative to the formalin fixation in sensitive Antigens.

Technical specifications:

WxDxH: 300x410x460mm

Weight: 31Kg

Power supply: 230V / 50Hz / 0.3KW

Article No. 038001