


Linear Stainer LS 72

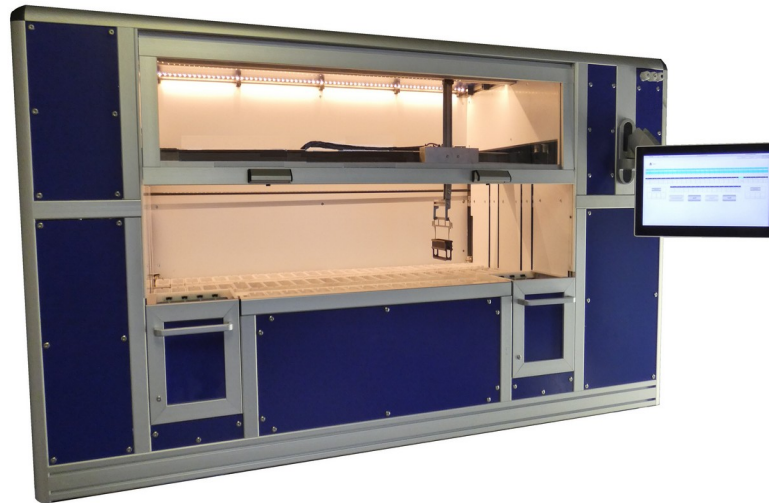

The growing need for space in histology laboratories, driven by the increasing number and size of devices, poses significant challenges for users. Our innovative solution addresses this issue by integrating two staining units within a single housing. Traditional staining machines typically feature only one arm for moving slide baskets, often limiting their capacity—handling more than five slide baskets simultaneously is uncommon due to varying staining protocols. As a result, users frequently require an additional staining unit, usually a linear stainer capable of managing multiple slide baskets at once. Our multifunctional staining machine seamlessly combines both device concepts. Designed to meet the highest standards, it offers a maximum throughput of 1,800 slides per hour.

Key Features:

- 26 positions in the linear staining row, with adjustable dwell times ranging from 45 seconds to 5 minutes
- Capacity to store up to 20 staining protocols
- Simultaneous staining with up to 40 slide baskets
- Loading and unloading stations, each accommodating 4 slide baskets
- Up to 6 freshwater cuvettes
- Adjustable agitation settings
- Optional built-in activated carbon filter or compatibility with existing exhaust air systems
- Adaptable to nearly all automatic cover-slippers
- Slide baskets available for 20 or 30 slides
- Integrated Reagent Management Program

Technical Data	Details
Dimensions (W x D x H)	200 x 75 x 120c (cm)
Connection values	230V / 50Hz / 1.2KW
Weight	230 kg